



Examination of blood coagulation


Homeostasis tests:



PT:



PTT:



Blood clotting or hemocoagulation, is the third stage of stopping bleeding. The first phase is formed first by reflex, later by humoral stimulation of **vasoconstriction**. n the second phase, there is **activation and adhesion of platelets** and the formation of a primary (unstable) platelet thrombus. The third phase then includes the activation of the plasma **coagulation cascade** and the formation of a fibrin network stabilizing the thrombus.

Indications for examination of blood coagulation

- **screening examination** – before invasive procedures (surgery, [[biopsy] etc.);
- suspected **bleeding or thrombophilic conditions**;
- **monitoring of anticoagulant treatment** – Warfarin;
- an auxiliary indicator in the **evaluation of hepatic proteosynthesis**.

Procedure

Blood is collected in a tube with EDTA or with citrate. These substances have the ability to bind calcium (form a complex with it) and thus prevent immediate blood clotting. The blood is sent without delay to the laboratory, where the fraction of blood cells is separated from the decalcified plasma by centrifugation. The required tests are then performed, e.g. aPTT, PT = Quick test, fibrinogen level or D-dimers.

APTT

APTT (*activated partial thromboplastin*, time) is a test of the internal and common pathway of hemocoagulation.

Execution

At a temperature of 37°C, **kaolin - kephalin complex** (kephalin - phospholipid from rabbit brain) is added to the decalcified plasma together with an excess of calcium ions. Here, kaolin represents a negatively charged surface that activates the internal pathway, and kephalin functions as a **partial thromboplastin** (it has nothing to do with tissue thromboplastin!), i.e. actually a substitute for **platelet phospholipid**, necessary for the activation of **factor X**. In this way, the internal and then the common pathway is activated. The time until the formation of a fibrin clot is measured.

Physiological values

- **25,9-40 s^[1]**
- APTR = value for comparison with standard plasma parameter is in the range of **0,83 - 1,3**
- APTT is mainly dependent on the initial events of hemocoagulation in the **internal pathway** (f. XII, XI, IX a VIII), because this part of coagulation takes place more slowly. Only in the second row does it capture the effectiveness of the factors of the common path, X, V, II and I.

Extension




- **hemophilia A** (f. VIII deficiency), **hemophilia B** (f. IX deficiency), **hemophilia C** (f. XI deficiency);
- **treatment with heparin** iv, but not with fractionated heparin sc - heparin increases the effectiveness of the physiological anticoagulant **antithrombin III**, which inhibits factors XII, XI, IX, X, II and VII. Except for f. VII, all these changes are included in the APTT. Heparin also directly affects factor X activity;
- **Warfarin treatment**;
- **von Willebrand's disease** – a hereditary disorder of vW factor, production, which is a carrier of factor VIII and thus affects its activity;
- **antiphospholipid syndrome** – presence of antibodies against phospholipids;;
- **consumptive coagulopathy** – DIC.

Quick test

right|náhled|Účinek warfarinu Quick's test, PT (prothrombin time, TT (thromboplastin time)) or thromboplastin time, is a test of the **external and common pathway**. Indicates the rate of conversion of prothrombin to thrombin due to the action of tissue thromboplastin (coagulation factor III).

Execution

Tissue thromboplastin is added to the **decalcified** plasma along with an excess of calcium ions. **In this way, the outer** and then the common path are activated. The time until the formation of the first fiber of the fibrin clot is measured. The results of the test are usually indicated by the **INR** index (*international normalized ratio*), i.e. the ratio of the patient's measured time and the normal value of the control plasma. Commercially available types of tissue thromboplastin (tissue factor III) differ among themselves. Therefore, the efficiency factor of the thromboplastin used is added to the equation for calculating the INR.

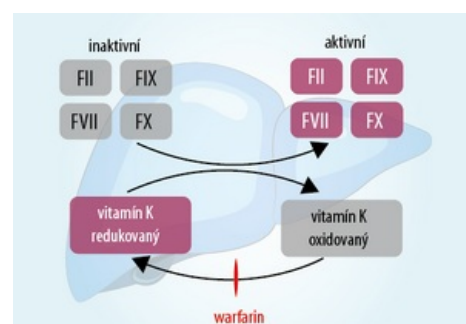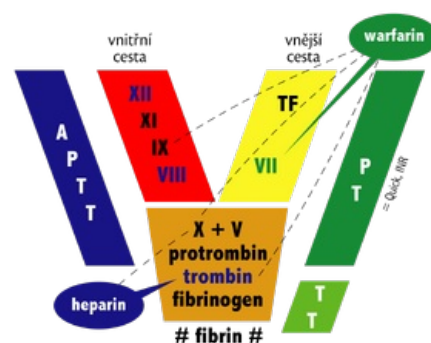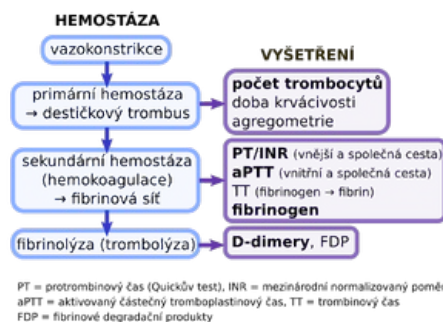

Physiological values

- **12-15 s**
- **INR: 0,8-1,2 (80-120 %)** ^[1]

The resulting time depends on the concentration of individual coagulation factors of the external and joint system. With increased blood coagulation, the INR is lower, with prolonged coagulation (e.g. during treatment with antivitamin K), the INR increases.

Extension

- **physiologically in newborns** – lack of factor VII;
- **Warfarin** therapy or other conditions with hypovitaminosis K;
- i.v. **heparin** therapy;
- severe **hepatic proteosynthesis disorder**;
- **consumptive coagulopathy** – DIC.

Further examination

The Lee-White method

- An indicative method of examining coagulation at the patient's bedside, which has minimal clinical use. The principle of the test consists in the spontaneous activation of the internal coagulation cascade *in vitro* (in a test tube). The result depends on **the state of the internal branch and the common** coagulation system.
- Normal time is: 5-15 min at room temperature.

Total fibrin degradation products

- FDP (fibrin degradation products) and **D-dimers** are specific for stabilized fibrin (norm below 160 µg/l).
- Immunological determination of fibrin degradation products follows acute phase proteins. It is an indicator of **activated fibrinolysis**, from a diagnostic point of view, primarily a marker of recently started coagulation with subsequent fibrinolysis. It is mainly used when thromboembolic disease is suspected. However, the high sensitivity of the method comes at the expense of its specificity. It can also be used to monitor treatment in which we are trying to induce clot dissolution in the patient.

Fibrinogen concentration

- It is referred to as **Fbg**. This is a determination of the plasma concentration of fibrinogen. The norm is between **1,5-4,5 g/l**.

Antitrombin III

- Determination of functional activity in plasma, normal range given by comparison with control plasma (70-100%). Primary or secondary (consumption in DIC) antithrombin III deficiency is a risk factor for thrombophilia.

Thrombin time

- Thrombin time (TT, *thrombin time*) is a direct examination of the **conversion of Fbg to fibrin**.
- **Execution:** thrombin is added to the decalcified plasma together with an excess of calcium. The time until the formation of a fibrin clot is measured.
- **The physiological range is 17-24 s.**
- Prolonged in case of dysfibrinogenemia, severe hypo- or afibrinogenemia, activated fibrinolysis and during treatment with heparin.

Euglobulin method for determination of fibrinolytic activity

- **The euglobulin fraction**, containing plasminogen, is precipitated from the plasma at the isoelectric point with acetic acid. The fractions are separated, dissolved and re-precipitated with calcium chloride. The plasmin present, freed from the effect of α₂-antiplasmin by precipitation, lyses the formed clot. When the fibrinolytic system is activated, more activated free plasmin is present and clot lysis is faster.
- **Physiological lysis time is 120-240 minutes.** Shorter times signal increased activation of the fibrinolytic system.

Activated coagulation time (ACT)

- Blood taken without anticoagulant is applied to a test tube with a contact activator (kaolin) and mixed in the device at 37 °C until the coagulation time is read. Routine for checking heparinization in extracorporeal circulation and hemodialysis. ACT of normal blood is approx. 150 s, with heparinization for long-term extracorporeal circulation or hemodialysis 180–300 s, with MO in cardiac surgery > 600 s.

Another examination

- examination of **protein C, protein S**, examination of **APC resistance** (genetically determined resistance of factor V to activated protein C, representing a risk factor for thromboembolic disease);
- examination of **antiphospholipid antibodies** – lupus anticoagulant, modified aPTT, anticardiolipin antibodies;
- examination of the **fibrinolytic system** – plasminogen activators (tPA, uPA, PAI-1), plasminogen;
- examination of the **activity of individual factors of the coagulation and fibrinolytic cascade** – dg. e.g. hemophilia A, B, C.

Evaluation of the examination of hemostasis

aPTT	PT (Quick)	Thrombocytes	Bleeding time	Common causes
				Vasculopathy, lack of f. XIII
			↑	Thrombocytopathy
		↓	↑	Thrombocytopenia
↑				Treatment with heparin , lack of f. VIII, IX, XI, XII
	↑			Lack of f. VII
↑	↑			Treatment with Warfarin and other hypovitaminosis K
↑			↑	von Willebrand disease
↑	↑	↓	↑	Impaired hepatic proteosynthesis, DIC (sepsis)

loosely based on E. Lechler

Links

Related Articles

WikiSkripta. *Krev* - *WikiSkripta* [online]. [cit. online]. <<https://www.wikiskripta.eu/w/Krev>>.

WikiSkripta. *Krevní obraz* - *WikiSkripta* [online]. [cit. online]. <https://www.wikiskripta.eu/w/Krevn%C3%AD_obraz>.

WikiSkripta. *Hemokoagulace* - *WikiSkripta* [online]. [cit. online]. <<https://www.wikiskripta.eu/w/Hemokoagulace>>.

WikiSkripta. *Vyšetření krvácivosti* - *WikiSkripta* [online]. [cit. online]. <https://www.wikiskripta.eu/w/Vy%C5%A1et%C5%99en%C3%AD_krv%C3%A1civosti>.

- Krev
- Krevní obraz
- Hemokoagulace
- Vyšetření krvácivosti

References

Maxdorf. *Velký lékařský slovník On-line: Normální laboratorní hodnoty dospělých* [online]. [cit. 2010-02-28]. <<https://lekarske.slovníky.cz/normalni-hodnoty>>.

1.

External links

- Fyziologický atlas – hemostáza (<http://www.physiome.cz/atlas/hemostaza/01/>)

Kategorie:Vložené články Kategorie:Hematologie Kategorie:Patofyziologie Kategorie:Fyziologie

